

Tubing Set CF


Tubing Set for CoolFlow™ pump for HF ablation catheter irrigation

The Tubing Set CF is being used together with the CoolFlow™ irrigation pump for irrigated HF-ablation catheters, to achieve a save an continuous cooling of the catheter tip. The Tubing Set CF consists of a drip chamber with spike to connect to an IV solution container, a short special tubing section to be placed within the irrigation pump, additional tubing to achieve

an overall length of 280 cm and a 3-way stopcock with Luer connections to connect to the ablation catheter. The special tubing section includes air bubble detection areas and fixation knobs. Please refer to the user manual of the CoolFlow™ irrigation pump, the HF generator and the HF ablation catheter for further user instructions.

Properties

- Compatible tubing set for Biosense Websters CoolFlow™ irrigation pump
- Drip chamber with spike to connect to an IV solution container
- 3-way stopcock with Luer connections to connect to ablation catheter
- Optimized for minimal wearout of the CoolFlow™ irrigation pump

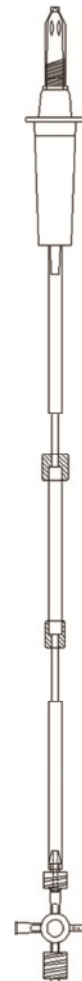

Order information

Tubing Set CF

Product	Irrigation pump	Catheter connection	IV container connection	Length	Order no.
Tubing Set CF #103-050	CoolFlow™ irrigation pump of Biosense Webster, Inc.	3-way stopcock with Luer connections	Drip chamber with spike	280 cm	39 1030 50

Order your Tubing Set CF

by phone: 0241-173518 | by fax: 0241-175627 | by e-mail: info@bispingmed.de

COOLFLOW™ is a registered trademark of Biosense Webster, Inc.